
waterpik™
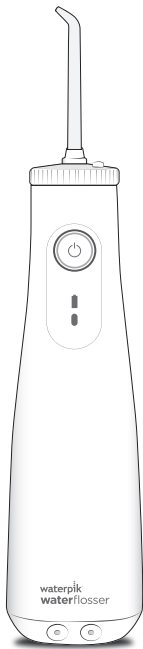
# Waterpik™ Water Flosser Model WF-10

Waterpik™ 沖牙機  
型號 WF-10

[www.waterpik.com](http://www.waterpik.com)

samcooper**design**

**DO NOT PRINT**

**DATE:** 02.07.2024

**CLIENT:** WATERPIK

**FILENAME:** 20033211-FAA\_WF-10W\_TraditionalL\_IM\_r6.indd

**DESCRIPTION:** WF-10 Instruction Manual

**LANGUAGES:** 2 [Eng, Trad Chinese]

**FILE TYPE:** Adobe InDesign CC

**SCD JOB# / PREPARED BY:** 23WPOH7979 / BS / BT

**CLIENT ROUND:** 6

**SCALE:** 100%

**UPC:** N/A

**SIZE:** 14.0 x 4.5 flat; 7.0 x 4.5 folded

**FORM NUMBER:** 20033211-F AA

**INK:** Process Black

**NOTES:**

# IMPORTANT SAFEGUARDS

## READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### DANGER:

#### To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place or drop the product into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not handle the user supplied USB power supply (power adapter) with wet hands.

- Unplug the user supplied USB power supply (power adapter) when not in use.
- Check the charging cord and user supplied USB power supply (power adapter) for damage before the first use and during the life of the product.

### WARNING:

#### To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only for its intended use as described in product literature or as recommended by your dental professional.

- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Water Pik, Inc. at 1-800-525-2774 or visit [www.waterpik.com](http://www.waterpik.com).
- A user supplied USB-A wall charger (power adapter), safety certified (IPX4, Class II) with a rated output of 5.0 VDC, 1.0 A is required to charge this product.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Keep the charging cord away from heated surfaces.
- Never drop or insert any object into any hose or opening.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Fill reservoir with warm water or other dental professional recommended solutions only.
- Only use tips and accessories recommended by Water Pik, Inc. as identified in this instruction manual.
- If the charging cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.

- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Do not direct water from the Water Flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- Do not use the Water Flosser for more than five minutes in each two-hour period.

- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
- There are no user serviceable electrical parts in this product and it does not require maintenance.
- Remove the charging cord from the unit prior to filling the reservoir.
- Ensure contacts on the cord, the unit and the power adapter are clean and dry while the unit is being charged.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet. Contact your doctor with any questions.
- Do not expose the battery or appliance to excessive temperatures.

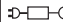
**SAVE THESE  
INSTRUCTIONS**

Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik™ Water Flosser	9
Battery Care, Battery Disposal/Recycling, Battery Removal	10
Limited One-Year Warranty	11


## SYMBOLS GLOSSARY

(Not all symbols apply to all products)


Explanation of symbols used on Water Pik, Inc. product labels and packaging. Applicable symbols for a specific product are found on that product's labels and packaging.

Symbol	Definition	Symbol	Definition	Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Detachable Power Supply		Do not use this device in a shower, bathtub and do not immerse in any liquid.		Regulatory Compliance Mark
	Manufacturer		Direct Current (DC) Voltage		Authorized Representative in the European Community		Separate collection for waste electrical and electronic equipment.
IPX7 or IPX4	Ingress Protection Code		Alternating Current (AC) Voltage		CE marking of conformity		Contains a recyclable Li-ion batter
	IEC Class II Equipment		Do not immerse this device in any liquid.		UK Conformity Assessed		Contains a recyclable Ni-MH battery

1. Tip rotator
2. Power button
3. Mode/pressure button
4. Battery charge indicator
5. Charging contacts
6. Tip
7. Tip eject button
8. Reservoir fill door
9. Reservoir

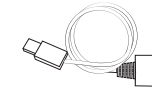
Charging Cable  
20026734-1 White
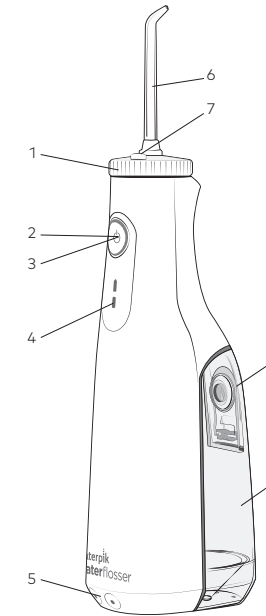
## TIPS




- |  |  |
|--|--|
| <b>A. Orthodontic Tip*</b><br>- Braces/General use                                   | <b>D. Classic Jet Tip**</b><br>- General use                                       |
| <b>B. Pik Pocket™ Tip*</b><br>- Periodontal Pockets/Furcations                       | <b>E. Tongue Cleaner**</b><br>- Fresher breath                                     |
| <b>C. Plaque Seeker™ Tip*</b><br>- Implants/Crowns/Bridges/<br>Retainers/General use | <b>F. Travel Plug</b><br>- Keeps excess water from<br>draining out while traveling |

## NOT ALL TIPS INCLUDED IN ALL MODELS

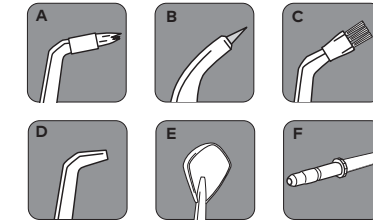\*Replace Every 3 months \*\*Replace Every 6 months

To find detailed information about tips, other accessories, and product usage please visit our website at [www.waterpik.com](http://www.waterpik.com) or email at [InternationalService@waterpik.com](mailto:InternationalService@waterpik.com).
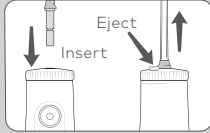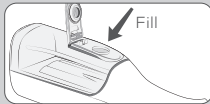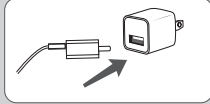
### Charging Unit Prior to First Use

Connect the USB charging cable to USB wall charger.

**NOTE:** See Important Safeguards for required USB wall charger specifications.

Plug the USB wall charger into a wall outlet and connect the magnetic charging cable to the unit and charge for 4 hours.

**Ensure contacts on BOTH the cord and unit are clean and dry before connecting the charger.**

The LED light will display a single, rapid blink every 3 seconds when unit is fully charged and the charger is connected. A full charge will last for approximately 1 week of typical use. When the batteries are in need of charging (approximately 30% charge left) the charge indicator light will blink 8 times quickly after use to indicate that you should recharge the unit. A complete recharge will take 3-4 hours. If you use your Water Flosser once a day or less, it is recommended you charge it once a week. If you use it twice a day or more, charge the unit more frequently.



### Filling the Reservoir

Remove the magnetic power cord from the unit prior to filling reservoir and product use. Lift flip top on reservoir and fill the reservoir with lukewarm water.

### Inserting and Removing Tips

Insert tip into the center of the knob at the top of the Water Flosser handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle.

**CAUTION: Do not press tip eject button while unit is running.**
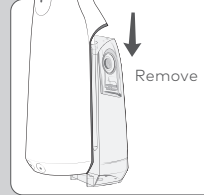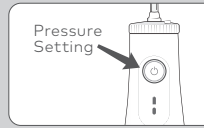
### Adjusting the Pressure Setting

- 2 pressure settings;
- Press power button to power on and switch settings.
- LED: Blink 3 times quickly every second:

Slow glow: the batteries are charging.

Single, rapid blink every 3 seconds: batteries are fully charged.

### Recommended Technique

Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON.


For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

### When Finished

Turn the unit off. Empty any liquid left in the reservoir.

If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove.

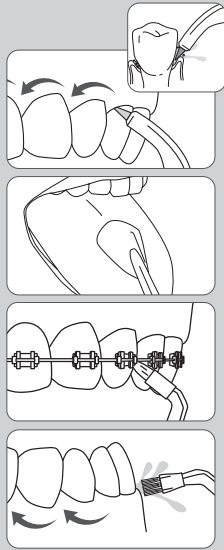**NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.**

#### TIP USAGE

##### Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

##### Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

##### Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

##### Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

#### Using Mouthwash and Other Solutions

Your Waterpik™ Water Flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions remove the reservoir (see "When Finished" section) and add from top opening on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

#### Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

It is recommended that the reservoir is removed and cleaned weekly.

#### Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.

### Battery Care

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage.

Do not expose the battery or appliance to excessive temperatures.

### Battery Disposal/Recycling




This product contains a **rechargeable Nickel Metal Hydride**. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations.

Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

Please contact your local waste administrative department for any questions.

**Note: The rechargeable battery is not replaceable, remove only at end of product life.** This process is not reversible and voids the product warranty.
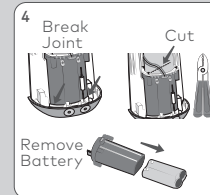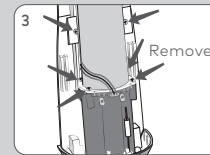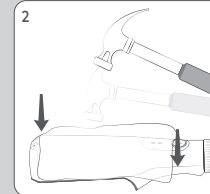
### Battery Removal

To remove the rechargeable battery for recycling, please discharge the battery fully before removing.

Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work.

1. Power off the device, remove reservoir.
2. Wrap in a towel and hit with a hammer to break the joint between the front, rear, and lower housings.
3. Remove 2 screws holding pump assembly to rear cover and 4 screws holding battery housing to pump assembly.
4. Carefully break the joint between the battery housing and rear cover.  
Cut battery wires (careful not to short positive and negative wires together) and remove battery.
5. Recycle/dispose of the battery at your nearest recycling center or take to a hazardous waste collection facility to comply with local and state regulations.

#### WARNING:

**DO NOT cut both battery wires at the same time.**

**DO NOT disassemble, puncture, hit, squeeze or throw the battery into fire.**

**Service Maintenance**

Waterpik™ Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all accessories go to [www.waterpik.com](http://www.waterpik.com) or email [InternationalService@waterpik.com](mailto:InternationalService@waterpik.com).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit. Replacement tips can be purchased from [www.waterpik.com](http://www.waterpik.com) or email [InternationalService@waterpik.com](mailto:InternationalService@waterpik.com).

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 1 year from the date of purchase.

**Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim.** Water Pik, Inc. will replace at its discretion the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips, etc. **Use with improper voltage will damage product and VOID warranty.**

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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©2024 Water Pik, Inc.
**Water Pik, Inc., a subsidiary of
Church & Dwight Co., Inc.**